



HANDY CARE BIDET

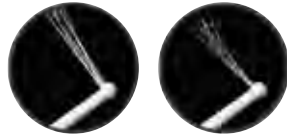

H: Strong Mode L: Soft Mode

Clean when & where you need.
Strong and Soft Mode, Swift Mode as you need

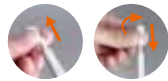

Waterpipe Port 180° adjustment for complete cleaning



1 minute automatically shut down, Energy Saving


Waterproof IPX6, easily clean and clear itself



Usage instructions:

Turn and draw bidet body from water reservoir; Fill with clean water (145ml), Fixed the bidet body with water reservoir
Open waterpipe in a proper way Starting the Strong or Soft mode after 3 seconds press power. Clean the bidet body after using.

Function

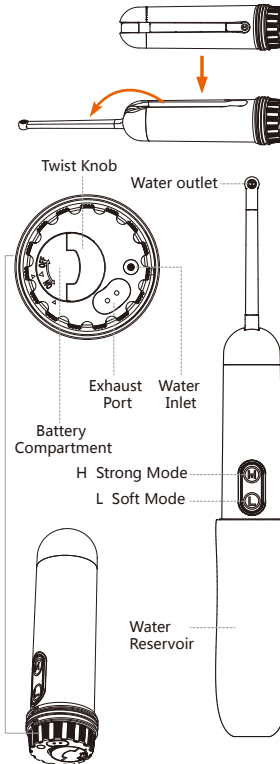

Technical Specifications:

Size	Φ47 X 156mm
Weight	215g
Reservoir Volume	145mL
Run Time(Soft Mode)	≥ 180 Times
Displacement(Soft Mode)	145mL/Minute
Run Time(Strong Mode)	≥ 120 Times
Displacement(Strong Mode)	300mL/Minute
Waterproof	IPX6
Noise	≤65dB
Working Output Current	< 200ma
Motor	DC 3.0V YYP370-3C2
Power Source	2 x AAA Battery
Working Mode	L Soft Mode-H Strong Mode

Cautions:

1. This bidet is intended for external use only.
2. Avoid filling with water above 40°C.
3. Clean the bidet body, reservoir and waterpipe on a regular basis to avoid harmful bacteria.
4. For proper cleaning use the bidet spray in a circular motion.
5. Prevent the device from falling, knocking and forced vibration.
6. Do not detach and change any internal parts.
7. Remove the battery if haven't used in a long time.
8. This product is not a toy. And contains small parts not suitable for children under 5 years of age